
neurotic Documentation

Release 1.0.0

Jeffrey Gill

Jul 27, 2019

Table of Contents


1	Overview	3
2	Installation	5
3	Getting Started	7
4	Configuring Metadata	9
4.1	Top-Level Organization	9
4.2	Specifying Data Locations	9
4.3	Remote Data Available for Download	10
4.4	Global Configuration Settings	11
4.5	Video Synchronization Parameters	12
4.6	Plot Parameters	13
4.7	Epoch Encoder Parameters	13
4.8	Filters	14
4.9	Amplitude Discriminators	14
4.10	tridesclous Spike Sorting Results	15
4.11	A Complete Example	16
5	API Reference Guide	19
5.1	neurotic.datasets.data	19
5.2	neurotic.datasets.download	20
5.3	neurotic.datasets.ftppath	20
5.4	neurotic.datasets.metadata	20
5.5	neurotic.gui.config	22
5.6	neurotic.gui.epochencoder	23
5.7	neurotic.gui.notebook	23
5.8	neurotic.gui.standalone	24
6	Release Notes	25
6.1	neurotic 1.0.0	25
6.2	neurotic 0.7.0	26
6.3	neurotic 0.6.0	26
6.4	neurotic 0.5.1	26
6.5	neurotic 0.5.0	26
6.6	neurotic 0.4.2	27
6.7	neurotic 0.4.1	27

6.8	neurotic 0.4.0	27
6.9	neurotic 0.3.0	27
6.10	neurotic 0.2.0	28
6.11	neurotic 0.1.1	28
6.12	neurotic 0.1.0	28
Python Module Index		29
Index		31

Curate, visualize, annotate, and share your behavioral ephys data using Python

Version: 1.0.0 (other versions)

neurotic is an app for Windows, macOS, and Linux that allows you to easily review and annotate your electrophysiology data and simultaneously captured video. It is an easy way to load your Neo-compatible data into ephyviewer without doing any programming. Share a single metadata file with your colleagues and they too will quickly be looking at the same datasets!

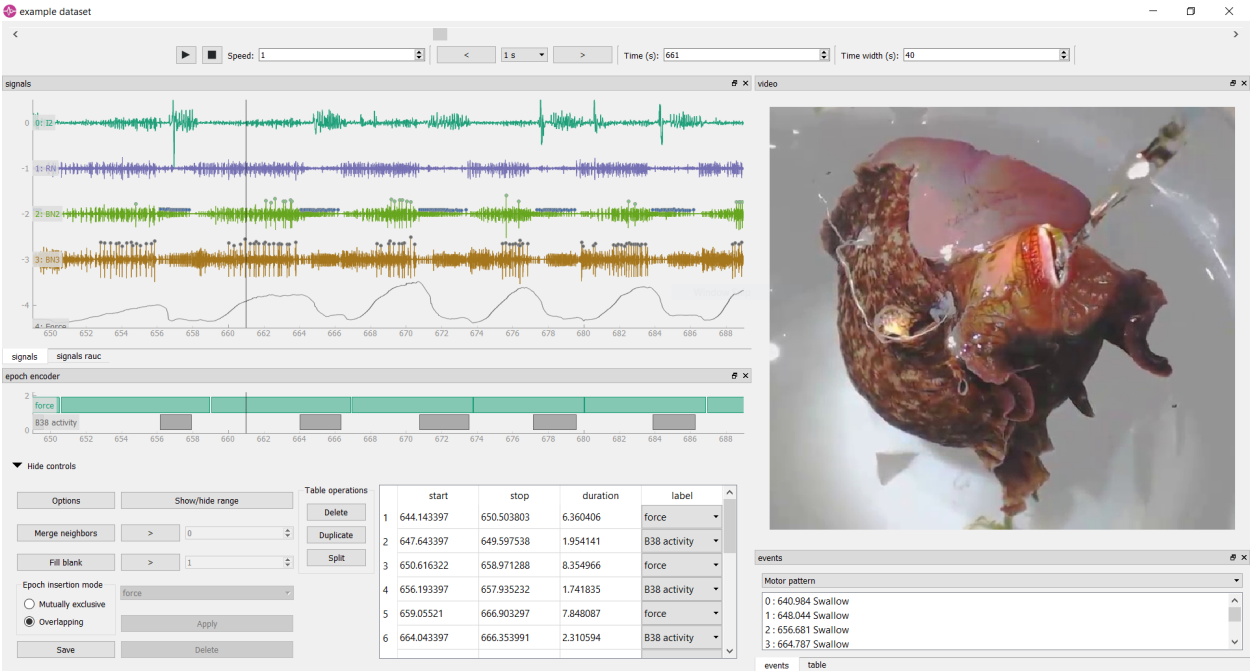

CHAPTER 1

Overview

To use **neurotic**, first organize your datasets in a *metadata file* like this (see *Configuring Metadata*):

```
my favorite dataset:
  description: This time it actually worked!

  data_dir:      C:\local_dir_containing_files
  remote_data_dir: http://myserver/remote_dir_containing_downloadable_files #_
↪ optional
  data_file:      data.axgx
  video_file:      video.mp4
  # etc

  video_offset: -3.4 # seconds between start of video and data acq
  epoch_encoder_possible_labels:
    - label01
    - label02
  plots:
    - channel: I2
      ylim: [-30, 30]
    - channel: RN
      ylim: [-60, 60]
    # etc

  filters: # used only if fast loading is off (lazy=False)
    - channel: Force
      lowpass: 50
    # etc


  amplitude_discriminators: # used only if fast loading is off (lazy=False)
    - name: B3 neuron
      channel: BN2
      units: uV
      amplitude: [50, 100]
    # etc
```

(continues on next page)

(continued from previous page)

```
another dataset:
# etc
```

Open your metadata file in **neurotic** and choose a dataset. If the data and video files aren't already on your local computer, the app can download them for you, even from a password-protected server. Finally, click launch and the app will use a standard viewer layout to display your data to you using **ephyviewer**.

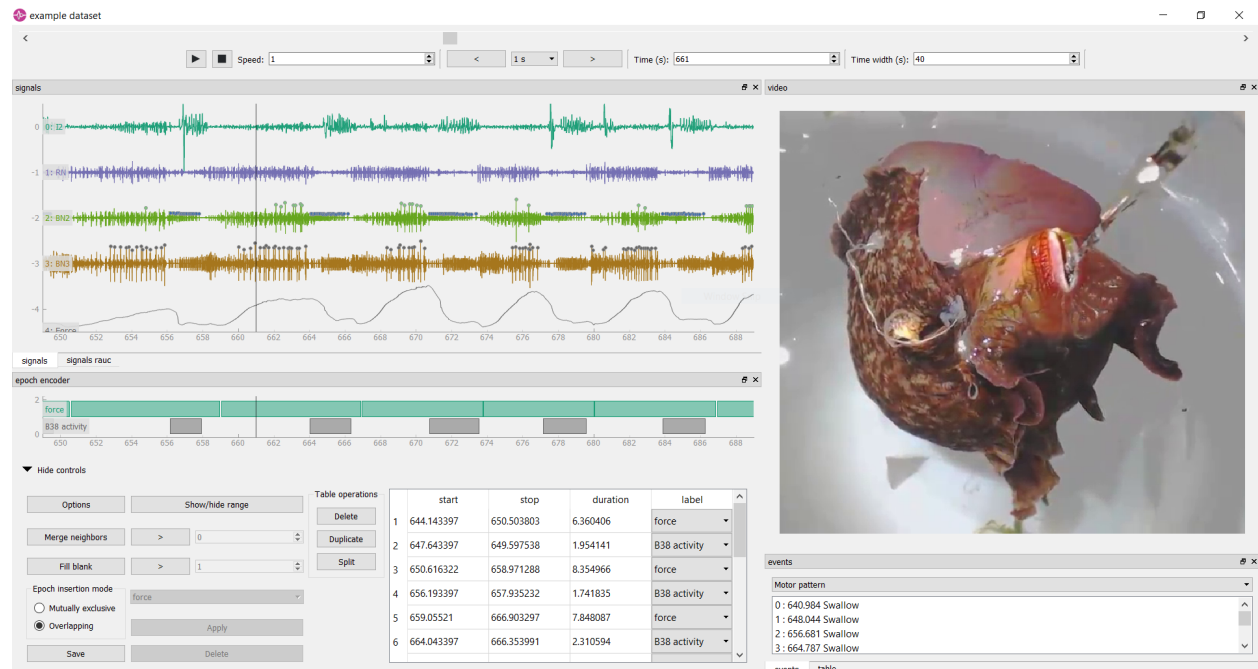

(Pictured above is a voracious *Aplysia californica* in the act of making the researcher very happy.)

The viewers are easy and intuitive to navigate (see [User Interface](#)):

- Pressing the play button will scroll through your data and video in real time, or at a higher or lower rate if the speed parameter is changed.
- The arrow/WASD keys allow you to step through time in variable increments.
- Jump to a time by clicking on an event in the event list or a table entry in the epoch encoder.
- To show more or less time at once, right-click and drag right or left to contract or expand time.
- Scroll the mouse wheel in the trace viewer or video viewer to zoom.
- The epoch encoder can be used to block out periods of time during which something interesting is happening for later review or further analysis (saved to a CSV file).
- All panels can be hidden, undocked, stacked, or repositioned on the fly.

Electrophysiologists will find this tool useful even if they don't need the video synchronization feature!

Portability is easy with neurotic! Use relative paths in your metadata file along with a remotely accessible data store such as [GIN](#) to make your metadata file fully portable. The same metadata file can be copied to a different computer, and downloaded files will automatically be saved to the right place. Data stores can be password protected and **neurotic** will prompt you for a user name and password. This makes it easy to share the **neurotic** experience with your colleagues!

CHAPTER 2

Installation

neurotic requires Python 3.6 or later. It also needs [PyAV](#), which is most easily installed using [conda](#):

```
conda install -c conda-forge av
```

PyAV must be installed manually, but all other dependencies will be installed with **neurotic**.

Install the latest release version of **neurotic** from [PyPI](#) using

```
pip install -U neurotic
```

or install the latest development version from [GitHub](#) using

```
pip install -U git+https://github.com/jpgill86/neurotic.git
```


CHAPTER 3

Getting Started

If you installed **neurotic** into a conda environment, first activate it:

```
conda activate <environment name>
```

Launch the app from the command line:

```
neurotic
```

A simple example is provided. Select the “example dataset”, download the associated data (~7 MB), and then click “Launch”. See [User Interface](#) for help with navigation.

Disabling “Fast loading” before launch will enable additional features including amplitude-threshold spike detection and signal filtering.

To inspect the metadata file associated with the example or to make changes to it, click “Edit metadata”. See [Configuring Metadata](#) for details about the format.

If you prefer Jupyter notebooks, you can launch an example notebook instead for experimenting with **neurotic**’s API:

```
neurotic --launch-example-notebook
```

The command line interface accepts other arguments too:

```
usage: neurotic [-h] [-V] [--no-lazy] [--thick-traces]
               [--theme {light,dark,original}] [--launch-example-notebook]
               [file] [dataset]

neurotic lets you curate, visualize, annotate, and share your behavioral ephys
data.

positional arguments:
  file                  the path to a metadata YAML file (default: an example
                        file)
  dataset               the name of a dataset in the metadata file to select
                        initially (default: the first entry in the metadata
```

(continues on next page)

(continued from previous page)

```
file)

optional arguments:
  -h, --help            show this help message and exit
  -V, --version          show program's version number and exit
  --no-lazy              do not use fast loading (default: use fast loading)
  --thick-traces         enable support for traces with thick lines, which has
                        a performance cost (default: disable thick line
                        support)
  --theme {light,dark,original}
                        a color theme for the GUI (default: light)
  --launch-example-notebook
                        launch Jupyter with an example notebook instead of
                        starting the standalone app (other args will be
                        ignored)
```

Configuring Metadata

To load your data with **neurotic**, you must organized them in one or more YAML files, called *metadata files*.

YAML files are very sensitive to punctuation and indentation, so mind those details carefully! Importantly, the tab character cannot be used for indentation; use spaces instead. There are [many free websites](#) that can validate YAML for you.

You may include comments in your metadata file, which should begin with #.

4.1 Top-Level Organization

Datasets listed within the same metadata file must be given unique names, which may include spaces. The special name `neurotic_config` is reserved for **neurotic** configuration settings and cannot be used for datasets.

In addition to names, a long description can be provided for each dataset.

Details pertaining to each dataset, including the description, are nested beneath the dataset name using indentation. You may need to use double quotes around names, descriptions, or other text if they contain special characters (such as : or #) or are composed only of numbers (such as a date).

```
experiment 2020-01-01:
  description: Both the name and description will be visible when neurotic loads_
↪the metadata
  # other details about this dataset will go here

my favorite dataset:
  description: This time it actually worked!
  # other details about this dataset will go here
```

4.2 Specifying Data Locations

Within a dataset's YAML block, paths to data and video files should be provided.

All files associated with a dataset should be collected into a single directory. A path to the local copy of this directory must be provided using the `data_dir` key. You may specify `data_dir` as an absolute path (e.g., `C:\Users\me\folder`) or as a path relative to the metadata file (e.g., `folder`).

Paths to individual files within the dataset are provided using keys listed below. These paths should be given relative to `data_dir`. If `data_dir` is flat (no subdirectories), these should be simply the file names. Only `data_file` is required.

Key	Description
<code>data_file</code>	A single Neo-compatible data file (required)
<code>video_file</code>	A video file that can be synchronized with <code>data_file</code>
<code>annotations_file</code>	A CSV file for read-only annotations
<code>epoch_encoder_file</code>	A CSV file for annotations writable by the epoch encoder
<code>tridesclous_file</code>	A CSV file output by tridesclous's <code>DataIO.export_spikes</code>

Note that the `annotations_file` must contain exactly 4 columns with these headers: “Start (s)”, “End (s)”, “Type”, and “Label”.

The `epoch_encoder_file` must contain exactly 3 columns with these headers: “Start (s)”, “End (s)”, and “Type”. (The fourth column is missing because ephyviewer’s epoch encoder is currently unable to attach notes to individual epochs; this may be improved upon in the future.)

The `tridesclous_file` is described in more detail in [tridesclous Spike Sorting Results](#).

4.3 Remote Data Available for Download

Data files must be stored on the local computer for **neurotic** to load them and display their contents. If the files are available for download from a remote server, **neurotic** can be configured to download them for you to the local directory specified by `data_dir` if the files aren’t there already.

Specify the URL to the directory containing the data on the remote server using `remote_data_dir`. **neurotic** expects the local `data_dir` and the `remote_data_dir` to have the same structure and will mirror the `remote_data_dir` in the local `data_dir` when you download data (not a complete mirror, just the specified files).

For an example, consider the following:

```
my favorite dataset:
  data_dir:          C:\Users\me\folder
  remote_data_dir:   http://myserver/remote_folder
  data_file:         data.axgx
  video_file:        video.mp4
```

With a metadata file like this, the file paths `data_file` and `video_file` are appended to `remote_data_dir` to obtain the complete URLs for downloading these files, and they will be saved to the local `data_dir`.

If you have many datasets hosted by the same server, you can specify the server URL just once using the special `remote_data_root` key, which should be nested under the reserved name `neurotic_config` outside of any dataset’s YAML block. This allows you to provide for each dataset a partial URL to a folder in `remote_data_dir` which is relative to `remote_data_root`. For example:

```
neurotic_config: # reserved name for global settings
  remote_data_root: http://myserver

my favorite dataset:
```

(continues on next page)

(continued from previous page)

```

data_dir:          C:\Users\me\folder1
remote_data_dir:   remote_folder1
data_file:         data.axgx
video_file:        video.mp4

another dataset:
data_dir:          C:\Users\me\folder2
remote_data_dir:   remote_folder2
data_file:         data.axgx
video_file:        video.mp4

```

Here, URLs to video files are composed by joining `remote_data_root + remote_data_dir + video_file`.

Recall that if `data_dir` is a relative path, it is assumed to be relative to the metadata file. In the example above, if the metadata file is located in `C:\Users\me`, the paths could be abbreviated:

```

neurotic_config:
  remote_data_root: http://myserver

my favorite dataset:
data_dir:          folder1
remote_data_dir:   remote_folder1
data_file:         data.axgx
video_file:        video.mp4

another dataset:
data_dir:          folder2
remote_data_dir:   remote_folder2
data_file:         data.axgx
video_file:        video.mp4

```

Note: Portability is easy with neurotic! Use relative paths in your metadata file along with a remotely accessible data store such as [GIN](#) to make your metadata file fully portable. The example above is a simple model of this style. A metadata file like this can be copied to a different computer, and downloaded files will automatically be saved to the right place. Data stores can be password protected and **neurotic** will prompt you for a user name and password. This makes it easy to share the **neurotic** experience with your colleagues!

4.4 Global Configuration Settings

The top-level name `neurotic_config` is reserved for configuration settings that apply to all datasets or to the app itself. Presently, only one configuration setting is implemented, but future versions of **neurotic** may add more under this name.

Key	Description
<code>remote_data_root</code>	A URL prepended to each <code>remote_data_dir</code> that is not already a full URL (i.e., does not already begin with a protocol scheme like <code>https://</code>)

For example:

```

neurotic_config:
  remote_data_root: http://myserver

```

(continues on next page)

```
my favorite dataset:
    # dataset details here
```

4.5 Video Synchronization Parameters

4.5.1 Constant Offset

If data acquisition began with some delay after video capture began, provide a negative value for `video_offset` equal to the delay in seconds. If video capture began after the start of data acquisition, use a positive value. A value of zero will have no effect.

neurotic warns users about the risk of async if `video_file` is given but `video_offset` is not. To eliminate this warning for videos that have no delay, provide zero.

4.5.2 Frame Rate Correction

If the average frame rate reported by the video file is a little fast or slow, you may notice your video and data going out of sync late in a long experiment. You can provide the `video_rate_correction` parameter to fix this. The reported average frame rate of the video file will be multiplied by this factor to obtain a new frame rate used for playback. A value less than 1 will decrease the frame rate and shift video events to later times. A value greater than 1 will increase the frame rate and shift video events to earlier times. A value of 1 has no effect.

You can obtain a good estimate of what value to use by taking the amount of time between two events in the video and dividing by the amount of time between the same two events according to the data record (seen, for example, as synchronization pulses or as movement artifacts).

4.5.3 Discrete Desynchronization Events

If you paused data acquisition during your experiment while video capture was continuous, you can use the `video_jumps` parameter to correct for these discrete desynchronization events, assuming you have some means of reconstructing the timing. For each pause, provide an ordered pair of numbers in seconds: The first is the time *according to data acquisition* (not according to the video) when the pause occurred, and the second is the duration of the pause during which the video kept rolling.

For example:

```
my favorite dataset:
    video_file: video.mp4
    # etc

    video_jumps:
        # a list of ordered pairs containing:
        # (1) time in seconds when paused occurred according to DAQ
        # (2) duration of pause in seconds
        - [60, 10]
        - [120, 10]
        - [240, 10]
```

These values could correct for three 10-second pauses occurring at times 1:00, 2:00, 3:00 according to the DAQ, which would correspond to times 1:00, 2:10, 3:20 according to the video. The extra video frames captured during the pauses will be excised from playback so that the data and video remain synced.

neurotic will automatically suggest values for `video_jumps` if it reads an AxoGraph file that contains stops and restarts (only if `video_jumps` is not already specified).

4.6 Plot Parameters

Use the `plots` parameter to specify which signal channels from `data_file` you want plotted and how to scale them.

Consider the following example, and notice the use of hyphens and indentation for each channel.

```
my favorite dataset:
  data_file: data.axgx
  # etc

  plots:
    - channel: Extracellular
      ylabel: Buccal nerve 2 (BN2)
      units: uV
      ylim: [-150, 150]

    - channel: Intracellular
      ylabel: B3 neuron
      units: mV
      ylim: [-100, 50]

    - channel: Force
      units: mN
      ylim: [-10, 500]
```

This would plot the “Extracellular”, “Intracellular”, and “Force” channels from the `data_file` in the given order. `ylabel` is used to relabel a channel and is optional. The `units` and `ylim` parameters are used together to scale each signal such that the given range fits neatly between the traces above and below it. If `units` is not given, they are assumed to be microvolts for voltage signals and millinewtons for force signals. If `ylim` is not given, they default to `[-120, 120]` for voltages and `[-10, 300]` for forces.

If `plots` is not provided, all channels are plotted using the default ranges, except for channels that match these patterns: “Analog Input #*” and “Clock”. Channels with these names can be plotted if given explicitly by `plots`.

The amount of time initially visible can be specified in seconds with `t_width`.

4.7 Epoch Encoder Parameters

The labels available to the epoch encoder must be specified ahead of time using `epoch_encoder_possible_labels` (this is a current limitation of ephyviewer that may eventually be improved upon).

For example:

```
my favorite dataset:
  epoch_encoder_file: epoch-encoder.csv
  # etc

  epoch_encoder_possible_labels:
    - label1
```

(continues on next page)

(continued from previous page)

- label12
- label13

4.8 Filters

Highpass, lowpass, and bandpass filtering can be applied to signals using the `filters` parameter. Note that filters are only applied if fast loading is off (`lazy=False`).

Consider the following example, and notice the use of hyphens and indentation for each filter.

```
my favorite dataset:
  data_file: data.axgx
  # etc

  filters: # used only if fast loading is off (lazy=False)

    - channel: Extracellular
      highpass: 300 # Hz
      lowpass: 500 # Hz

    - channel: Intracellular
      highpass: 300 # Hz

    - channel: Force
      lowpass: 50 # Hz
```

Filter cutoffs are given in hertz. Combining `highpass` and `lowpass` provides bandpass filtering.

4.9 Amplitude Discriminators

Spikes with peaks that fall within amplitude windows given by `amplitude_discriminators` can be automatically detected by **neurotic** on the basis of amplitude alone. Note that amplitude discriminators are only applied if fast loading is off (`lazy=False`).

Detected spikes are indicated on the signals with markers, and spike trains are displayed in a raster plot.

In addition to restricting spike detection for a given unit to an amplitude window, detection can also be limited in time to overlap with epochs with a given label.

Consider the following example, and notice the use of hyphens and indentation for each amplitude discriminator.

```
my favorite dataset:
  data_file: data.axgx
  # etc

  amplitude_discriminators: # used only if fast loading is off (lazy=False)

    - name: Unit 1
      channel: Extracellular
      units: uV
      amplitude: [50, 150]

    - name: Unit 2
```

(continues on next page)

(continued from previous page)

```

channel: Extracellular
units: uV
amplitude: [20, 50]
epoch: Unit 2 activity

```

Here two units are detected on the same channel with different amplitude windows. Any peaks between 50 and 150 microvolts on the “Extracellular” channel will be tagged as a spike belonging to “Unit 1”. The discriminator for “Unit 2” provides the optional `epoch` parameter. This restricts detection of “Unit 2” to spikes within the amplitude window that occur at the same time as epochs labeled “Unit 2 activity”. These epochs can be created by the epoch encoder (reload required to rerun spike detection at launch-time), specified in the read-only `annotations_file`, or even be contained in the `data_file` if the format supports epochs.

Amplitude windows are permitted to be negative.

4.10 tridesclous Spike Sorting Results

`tridesclous` is a sophisticated spike sorting toolkit. The results of a sorting process can be exported to a CSV file using `tridesclous's DataIO.export_spikes` function. This file contains two columns: the first is the sample index of a spike, and the second is the ID for a cluster of spikes. If this file is specified with `tridesclous_file`, then a mapping from the cluster IDs to channels must be provided with `tridesclous_channels`.

In the following example, notice the lack of hyphens:

```

my favorite dataset:
  data_file: data.axgx
  tridesclous_file: spikes.csv
  # etc

  tridesclous_channels:
    0: [Channel A, Channel B]
    1: [Channel A]
    2: [Channel B]
    3: [Channel B]
    # etc

```

Here numeric cluster IDs are paired with a list of channels found in `data_file` on which the spikes were detected.

To show only a subset of clusters or to merge clusters, add the `tridesclous_merge` parameter.

In this example, note again the punctuation:

```

my favorite dataset:
  data_file: data.axgx
  tridesclous_file: spikes.csv
  # etc

  tridesclous_channels:
    0: [Channel A, Channel B]
    1: [Channel A]
    2: [Channel B]
    3: [Channel B]
    # etc

  tridesclous_merge:
    - [0, 1]
    - [3]

```

Now clusters 0 and 1 are combined into a single unit, and only that unit and cluster 3 are plotted; cluster 2 has been discarded.

4.11 A Complete Example

These are the contents of the example metadata file that ships with **neurotic**, which can be loaded by running **neurotic** from the command line without arguments:

```
example dataset:
description: This is an example data set

# these data are a subset of Jeffrey Gill's dataset 2018-06-21_IN-VIVO_JG-08 002
data_dir:      example-data
remote_data_dir: https://web.gin.g-node.org/jpgill86/neurotic-data/raw/master/
→examples/example-data
data_file:      data.axgx
video_file:     video.mp4
annotations_file: annotations.csv
epoch_encoder_file: epoch-encoder.csv

video_offset: 640.3 # seconds

epoch_encoder_possible_labels:
- force
- B38 activity

plots:
- channel: I2
  units: uV
  ylim: [-30, 30]

- channel: RN
  units: uV
  ylim: [-60, 60]

- channel: BN2
  units: uV
  ylim: [-120, 120]

- channel: BN3
  units: uV
  ylim: [-150, 150]

- channel: Force
  units: mN
  ylim: [-10, 300]

filters: # used only if fast loading is off (lazy=False)

- channel: I2
  lowpass: 100 # Hz

- channel: Force
  lowpass: 50 # Hz

amplitude_discriminators: # used only if fast loading is off (lazy=False)
```

(continues on next page)

(continued from previous page)

```
- name: B3
  channel: BN2
  units: uV
  amplitude: [50, 150]

- name: B38
  channel: BN2
  units: uV
  amplitude: [17, 26]
  epoch: B38 activity

- name: B4/B5
  channel: BN3
  units: uV
  amplitude: [85, 200]
```

API Reference Guide

In addition to using **neurotic** as a standalone app, you can also leverage its API in your own code.

The core of the API consists of two classes and one function:

- `neurotic.datasets.metadata.MetadataSelector`: Read metadata files, download datasets
- `neurotic.datasets.data.load_dataset()`: Read datasets, apply filters and spike detection
- `neurotic.gui.config.EphyviewerConfigurator`: Launch ephyviewer

All public package contents are automatically imported directly into the `neurotic` namespace. This means that a class like `neurotic.datasets.metadata.MetadataSelector` can be accessed more compactly as `neurotic.MetadataSelector`.

5.1 `neurotic.datasets.data`

The `neurotic.datasets.data` module implements a function for loading a dataset from selected metadata.

`neurotic.datasets.data.load_dataset(metadata, lazy=False, signal_group_mode='split-all', filter_events_from_epochs=False)`

Load a dataset.

`metadata` may be a `MetadataSelector` or a simple dictionary containing the appropriate data.

The `data_file` in `metadata` is read into a `NeoBlock` using an automatically detected `neo.io` class if `lazy=False` or a `neo.rawio` class if `lazy=True`.

Epochs and events loaded from `annotations_file` and `epoch_encoder_file` and spike trains loaded from `tridesclous_file` are added to the `NeoBlock`.

If `lazy=False`, filters given in `metadata` are applied to the signals and amplitude discriminators are run to detect spikes.

5.2 neurotic.datasets.download

The `neurotic.datasets.download` module implements a general purpose download function that handles connecting to remote servers, performing authentication, and downloading files with progress reporting. The function handles various errors and will automatically reprompt the user for login credentials if a bad user name or password is given.

The module installs an `urllib.request.HTTPBasicAuthHandler` and a `neurotic.datasets.ftpauth.FTPBasicAuthHandler` at import time.

```
neurotic.datasets.download.download(url,          local_file,          overwrite_existing=False,
                                     show_progress=True, bytes_per_chunk=8192)
```

Download a file.

5.3 neurotic.datasets.ftpauth

The `neurotic.datasets.ftpauth` module implements a `urllib.request`-compatible FTP handler that prompts for and remembers passwords.

class `neurotic.datasets.ftpauth.FTPBasicAuthHandler` (*password_mgr=None*)

This subclass of `urllib.request.FTPHandler` implements basic authentication management for FTP connections. Like `urllib.request.HTTPBasicAuthHandler`, this handler for FTP connections has a password manager that it checks for login credentials before connecting to a server.

This subclass also ensures that file size is included in the response header, which can fail for some FTP servers if the original `FTPHandler` is used.

This handler can be installed globally in a Python session so that calls to `urllib.request.urlopen('ftp://...')` will use it automatically:

```
>>> handler = FTPBasicAuthHandler()
>>> handler.add_password(None, uri, user, passwd)  # realm must be None
>>> opener = urllib.request.build_opener(handler)
>>> urllib.request.install_opener(opener)
```

`neurotic.datasets.ftpauth.setup_ftpauth()`

Install `neurotic.datasets.ftpauth.FTPBasicAuthHandler` as the global default FTP handler.

Note that `urllib.request.install_opener()` used here will remove all other non-default handlers installed in a different opener, such as an `urllib.request.HTTPBasicAuthHandler`.

5.4 neurotic.datasets.metadata

The `neurotic.datasets.metadata` module implements a class for reading metadata files.

class `neurotic.datasets.metadata.MetadataSelector` (*file=None, local_data_root=None, remote_data_root=None, initial_selection=None*)

A class for managing metadata.

A metadata file can be specified at initialization, in which case it is read immediately. The file contents are stored as a dictionary in `all_metadata`.

```
>>> metadata = MetadataSelector(file='metadata.yml')
>>> print(metadata.all_metadata)
```


File contents can be reloaded after they have been changed, or after changing file, using the `load()` method.

```
>>> metadata = MetadataSelector()
>>> metadata.file = 'metadata.yml'
>>> metadata.load()
```

A particular metadata set contained within the file can be selected at initialization with `initial_selection` or later using the `select()` method. After making a selection, the selected metadata set is accessible at `metadata.selected_metadata`, e.g.

```
>>> metadata = MetadataSelector(file='metadata.yml')
>>> metadata.select('Data Set 5')
>>> print(metadata.selected_metadata['data_file'])
```

A compact indexing method is implemented that allows the selected metadata set to be accessed directly, e.g.

```
>>> print(metadata['data_file'])
```

This allows the `MetadataSelector` to be passed to functions expecting a simple dictionary corresponding to a single metadata set, and the selected metadata set will be used automatically.

Files associated with the selected metadata set can be downloaded individually or all together, e.g.

```
>>> metadata.download('video_file')
```

or

```
>>> metadata.download_all_data_files()
```

The absolute path to a local file or the full URL to a remote file associated with the selected metadata set can be resolved with the `abs_path()` and `abs_url()` methods, e.g.

```
>>> print(metadata.abs_path('data_file'))
>>> print(metadata.abs_url('data_file'))
```

abs_path(file)

Convert the relative path of `file` to an absolute path using `data_dir`.

abs_url(file)

Convert the relative path of `file` to a full URL using `remote_data_dir`.

all_metadata = None

A dictionary containing the entire file contents, set by `load()`.

download(file, **kwargs)

Download a file associated with the selected metadata set.

See `neurotic.datasets.download.download()` for possible keyword arguments.

download_all_data_files(kwargs)**

Download all files associated with the selected metadata set.

See `neurotic.datasets.download.download()` for possible keyword arguments.

load()

Read the metadata file.

select(selection)

Select a metadata set.

selected_metadata

The access point for the selected metadata set.

5.5 neurotic.gui.config

The *neurotic.gui.config* module implements a class for configuring and launching ephyviewer for a loaded dataset.

class neurotic.gui.config.EphyviewerConfigurator(*metadata, blk, rauc_sigs=None, lazy=False*)

A class for launching ephyviewer for a dataset with configurable viewers.

At initialization, invalid viewers are automatically disabled (e.g., the video viewer is disabled if `video_file` is not given in `metadata`). Viewers can be hidden or shown before launch using the built-in methods. Valid viewer names are:

- `traces`
- `traces_rauc`
- `freqs`
- `spike_trains`
- `epochs`
- `epoch_encoder`
- `video`
- `event_list`
- `data_frame`

launch_ephyviewer() is provided for starting a new Qt app and launching the ephyviewer main window all at once. *create_ephyviewer_window()* generates just the ephyviewer window and should be used if there is already a Qt app running.

create_ephyviewer_window(*theme='light', support_increased_line_width=False*)

Load data into each ephyviewer viewer and return the main window.

disable(*name*)

Disable the viewer name.

enable(*name*)

Enable the viewer name.

hide(*name*)

Hide the viewer name.

hide_all()

Hide all viewers.

is_enabled(*name*)

Return whether the viewer name is enabled.

is_shown(*name*)

Return whether the viewer name is shown.

launch_ephyviewer(*theme='light', support_increased_line_width=False*)

Start a Qt app and create an ephyviewer window.

show (*name*)
 Show the viewer name.

show_all ()
 Show all viewers.

5.6 neurotic.gui.epochencoder

The `neurotic.gui.epochencoder` module implements a subclass of `ephyviewer.datasources.epochs.WritableEpochSource`.

```
class neurotic.gui.epochencoder.NeuroticWritableEpochSource (filename,           pos-
                                                             sible_labels,
                                                             color_labels=None,
                                                             channel_name="",
                                                             backup=True)
```

A subclass of `ephyviewer.datasources.epochs.WritableEpochSource` for custom CSV column formatting and automatic file backup.

5.7 neurotic.gui.notebook

The `neurotic.gui.notebook` module implements Jupyter notebook widget counterparts for the `MetadataSelector` and the `EphyviewerConfigurator`.

```
class neurotic.gui.notebook.MetadataSelectorWidget (file=None,           lo-
                                                    cal_data_root=None,           re-
                                                    mote_data_root=None,           ini-
                                                    tial_selection=None)
```

Interactive list box for Jupyter notebooks that allows the user to select which metadata set they would like to work with.

```
>>> metadata = MetadataSelectorWidget (file='metadata.yml')
>>> display(metadata)
```

After clicking on an item in the list, the selected metadata set is accessible at `metadata.selected_metadata`, e.g.

```
>>> metadata.selected_metadata['data_file']
```

A compact indexing method is implemented that allows the selected metadata set to be accessed directly, e.g.

```
>>> metadata['data_file']
```

This allows the `MetadataSelectorWidget` to be passed to functions expecting a simple dictionary corresponding to a single metadata set, and the selected metadata set will be used automatically.

```
class neurotic.gui.notebook.EphyviewerConfiguratorWidget (metadata,           blk,
                                                            rauc_sigs=None,
                                                            lazy=False)
```

Interactive button grid for Jupyter notebooks that allows the user to select which ephyviewer viewers they would like to display and then launch ephyviewer.

disable (*name*)
 Disable the viewer name.

enable (*name*)
Enable the viewer *name*.

hide (*name*)
Hide the viewer *name*.

show (*name*)
Show the viewer *name*.

5.8 neurotic.gui.standalone

The *neurotic.gui.standalone* module implements the main window of the app.

```
class neurotic.gui.standalone.MainWindow(file=None, initial_selection=None,  
                                          lazy=True, theme='light', sup-  
                                          port_increased_line_width=False)
```

The main window of the app.

6.1 neurotic 1.0.0

2019-07-27

First stable release!

6.1.1 Improvements

- Major API changes ([#104](#), [#100](#), [#106](#))
 - In preparation for this stable release, many formerly public classes and functions were made private. This was done to minimize the number of public classes/functions, which beginning with this release will be treated as stable APIs that are ideally modified only in backwards compatible ways. Users should trust that public classes and functions will not be changed without good reason and a major version bump.
- Many improvements to the documentation, including the addition of an [API Reference Guide](#)
- Add example Jupyter notebook and command line argument for launching it ([#108](#))
- Add file overwrite option to download functions ([#106](#))
- Reserve the metadata keyword `neurotic_config` for global parameters ([#93](#))
 - The `remote_data_root` key must now be nested under `neurotic_config`.

6.1.2 Bug fixes

- Fix crash when epoch encoder file contains labels not listed in metadata ([#97](#))
- Allow amplitude discriminators to be specified with arbitrary units ([#99](#))

6.2 neurotic 0.7.0

2019-07-21

6.2.1 Improvements

- New documentation hosted at Read the Docs: <https://neurotic.readthedocs.io>
- Add menu action for opening metadata in editor (#83)
- Add menu action for opening the selected data directory (#84)
- Add list of installed versions of dependencies and doc links to About window (#44, #65)

6.2.2 Bug fixes

- Fix files remaining locked after closing a fast-loaded window (#69)
- Fix launching from command line with bad metadata argument (#82)

6.3 neurotic 0.6.0

2019-07-10

6.3.1 Improvements

- Add a basic “About neurotic” window with version and website information (#38)
- Update logo (#39)
- Add keywords and project URLs to package metadata (#40)

6.4 neurotic 0.5.1

2019-07-09

6.4.1 Compatibility updates

- Compatibility update for RawIOs with non-zero offset (#37)

6.5 neurotic 0.5.0

2019-07-06

6.5.1 Improvements

- Support fast (lazy) loading in Neo < 0.8.0 (#35)
- Add “git.” and conditionally “.dirty” to dev local version identifier (#34)

6.6 neurotic 0.4.2

2019-07-06

6.6.1 Bug fixes

- Fix for EstimateVideoJumpTimes regression introduced in 0.4.0 (#33)

6.7 neurotic 0.4.1

2019-07-02

6.7.1 Compatibility updates

- Change sources of development versions of dependencies (#32)
- Compatibility update for scaling of raw signals (#31)

6.8 neurotic 0.4.0

2019-07-01

6.8.1 Improvements

- Show epochs imported from CSV files with zero duration in epoch viewer (#27)
- Show epochs/events imported from data file in events list/epoch viewer (#28)
- Alphabetize epoch and event channels by name (#29)

6.9 neurotic 0.3.0

2019-06-29

6.9.1 Improvements

- Remove dependency on ipywidgets by making notebook widgets optional (#25)
 - Notebook widget classes renamed: `MetadataSelector` → `MetadataSelectorWidget`,
`EphyviewerConfigurator` → `EphyviewerConfiguratorWidget`
- Add app description and screenshot to README (#22)
- Promote to beta status (#23)

6.10 neurotic 0.2.0

2019-06-28

6.10.1 Improvements

- Add basic command line arguments (#14)
- Add continuous integration with Travis CI for automated testing (#13)
- Add some tests (#15, #16)
- Migrate example data to GIN (#18)

6.10.2 Bug fixes

- Fix crash when downloading from a server that does not report file size (#17)
- Raise an exception if a Neo RawIO cannot be found for the data file (#12)

6.11 neurotic 0.1.1

2019-06-22

6.11.1 Bug fixes

- Fix various downloader errors (#7)

6.12 neurotic 0.1.0

2019-06-22

- First release

n

- `neurotic.datasets.data`, [19](#)
- `neurotic.datasets.download`, [20](#)
- `neurotic.datasets.ftpauth`, [20](#)
- `neurotic.datasets.metadata`, [20](#)
- `neurotic.gui.config`, [22](#)
- `neurotic.gui.epochencoder`, [23](#)
- `neurotic.gui.notebook`, [23](#)
- `neurotic.gui.standalone`, [24](#)

A

`abs_path()` (*neurotic.datasets.metadata.MetadataSelector*
method), 21

`abs_url()` (*neurotic.datasets.metadata.MetadataSelector*
method), 21

`all_metadata` (*neurotic.datasets.metadata.MetadataSelector*
attribute), 21

C

`create_ephyviewer_window()` (*neurotic.gui.config.EphyviewerConfigurator*
method), 22

D

`disable()` (*neurotic.gui.config.EphyviewerConfigurator*
method), 22

`disable()` (*neurotic.gui.notebook.EphyviewerConfiguratorWidget*
method), 23

`download()` (in module *neurotic.datasets.download*), 20

`download()` (*neurotic.datasets.metadata.MetadataSelector*
method), 21

`download_all_data_files()` (*neurotic.datasets.metadata.MetadataSelector*
method), 21

E

`enable()` (*neurotic.gui.config.EphyviewerConfigurator*
method), 22

`enable()` (*neurotic.gui.notebook.EphyviewerConfiguratorWidget*
method), 23

`EphyviewerConfigurator` (class in *neurotic.gui.config*), 22

`EphyviewerConfiguratorWidget` (class in *neurotic.gui.notebook*), 23

F

`FTPBasicAuthHandler` (class in *neurotic.datasets.ftpauth*), 20

H

`hide()` (*neurotic.gui.config.EphyviewerConfigurator*
method), 22

`hide()` (*neurotic.gui.notebook.EphyviewerConfiguratorWidget*
method), 24

`hide_all()` (*neurotic.gui.config.EphyviewerConfigurator*
method), 22

I

`is_enabled()` (*neurotic.gui.config.EphyviewerConfigurator*
method), 22

`is_shown()` (*neurotic.gui.config.EphyviewerConfigurator*
method), 22

L

`launch_ephyviewer()` (*neurotic.gui.config.EphyviewerConfigurator*
method), 22

`load()` (*neurotic.datasets.metadata.MetadataSelector*
method), 21

`load_dataset()` (in module *neurotic.datasets.data*), 19

M

`MainWindow` (class in *neurotic.gui.standalone*), 24

`MetadataSelector` (class in *neurotic.datasets.metadata*), 20

`MetadataSelectorWidget` (class in *neurotic.gui.notebook*), 23

N

`neurotic.datasets.data` (module), 19

`neurotic.datasets.download` (module), 20

`neurotic.datasets.ftpauth` (module), 20

`neurotic.datasets.metadata` (module), 20

`neurotic.gui.config` (module), 22

`neurotic.gui.epochencoder` (module), 23

`neurotic.gui.notebook` (module), 23

`neurotic.gui.standalone` (*module*), [24](#)

`NeuroticWritableEpochSource` (*class in neurotic.gui.epochencoder*), [23](#)

S

`select()` (*neurotic.datasets.metadata.MetadataSelector method*), [21](#)

`selected_metadata` (*neurotic.datasets.metadata.MetadataSelector attribute*), [21](#)

`setup_ftpauth()` (*in module neurotic.datasets.ftpauth*), [20](#)

`show()` (*neurotic.gui.config.EphyviewerConfigurator method*), [22](#)

`show()` (*neurotic.gui.notebook.EphyviewerConfiguratorWidget method*), [24](#)

`show_all()` (*neurotic.gui.config.EphyviewerConfigurator method*), [23](#)